
# FeelYou2

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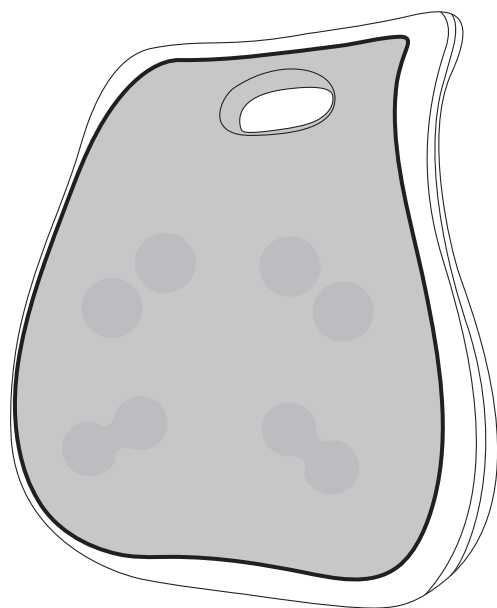Адресс производителя: D Building No. 100, JINFU ROAD,  
Chengnan Industrial Zone, TONG'AN District, Xiamen, China

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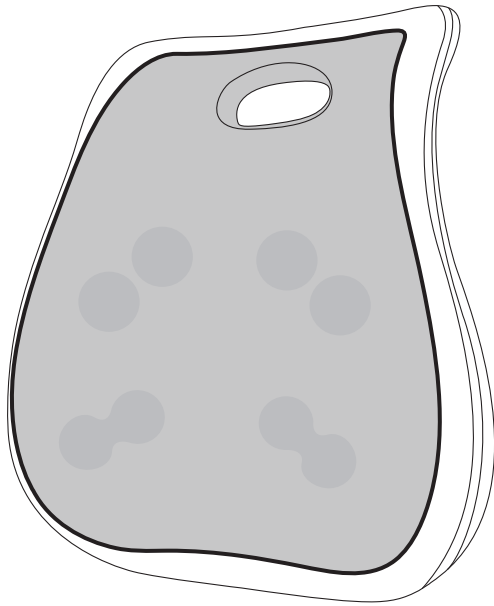
**User Manual**

# Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

**DANGER** - To reduce the risk of electric shock: Always unplug this product from the electrical outlet immediately after use and before cleaning

**WARNING** - To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not use and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Never operate the product with the air openings blocked. Keep the air openings away from lint, hair, and something like that.
- 4) Do not drop or insert any object into any opening.
- 5) Do not use outdoors.
- 6) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 7) Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
- 8) Do not carry this product by its supply cord or use cord as a handle.
- 9) Keep the cord away from heated surfaces. If the fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. please do not repair it by yourself, please contact the seller or manufacturer in time, avoid accidents or injuries.
- 10) Remove the power plug from the socket immediately after the power failure. Otherwise, it may cause electric shock, short circuit and fire.
- 11) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examinatio.
- 12) Special attention to components most susceptible to wear an advice to replace defective components immediately and/or keep the equipment out of use until repair.
- 13) Use this appliance only for its intended use as described in this manual.
- 14) The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- 15) The appliance is only to be used with the power supply unit provided with the appliance.
- 16) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

※ Children should be supervised to ensure that they do not play with the appliance.

※ Cleaning and user maintenance shall not be made by children without supervision.

17) This appliance contains batteries that are non-replaceable.

18) This appliance contains batteries that are only replaceable by skilled persons.

19) For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.

\* If you feel abnormal during use, please stop using immediately and consult your doctor. Otherwise may it be harmful to health.

1. People who are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implantable medical electronic devices.
2. People who are supposed to lose their sense of temperature.
3. People with malignant tumors.
4. People with heart disease.
5. A person who is restless early in pregnancy or shortly after giving birth.
6. One who has an abnormal or curved spine.
7. A person with a disturbance of perception caused by a high degree of peripheral circulation in diabetes.
8. Someone with a scar on their skin.
9. People who need quiet and people who are in poor health.
10. People with osteoporosis, severe fractures, sprain, and muscle pain, etc., acute (head) pain.
11. People with 38 °C temperature (fever).
12. (acute inflammatory symptoms, tiredness, chills, changes in blood pressure, etc.) a time of intense, debilitating stress.
13. People with tendonitis.
14. Others treated at medical facilities.

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using this product.

Only use indoors!

Do not use the device in damp rooms (e.g. when bathing or showering)

• Do not wash!

• Do not dry clean!

# SAVE THESE INSTRUCTIONS

## Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

## Introduction



Thank you for purchasing the Portable Back Massage Cushion for Lower Back. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.

### Package Contents

- 1 x Shiatsu Massage Cushion
- 1 x User Manual
- 1 x Car Charger
- 1 x Charger

### Product Diagram

1. Handle
2. Massage Node
3. Zipper
4. Power Button
5. DC Connector
6. Power Plug
7. Cable
8. Appliance Plug
9. Cigarette Lighter
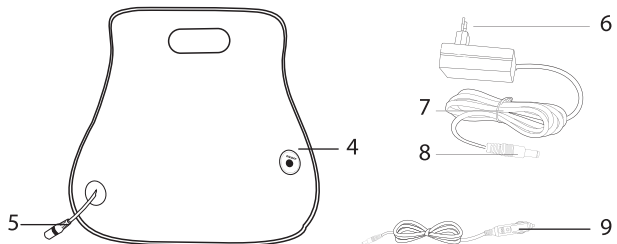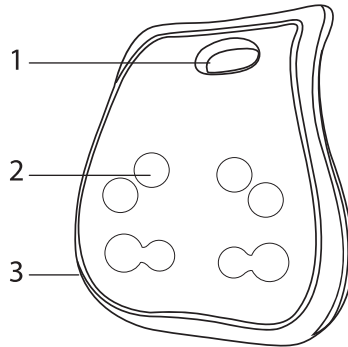
# How to Operate

1. Plug the power adapter into the 220V / 50Hz outlet.
2. Press the power button to active the massage and heating function. Press it again to turn off the heating function. The selection of heating is cycled.
3. Long press the power button for 3 seconds to turn off the device.
4. The device will change the massage direction automatically per 1 min. (Outward / Inward)
5. The device will stop automatically after 20 mins.
6. Please turn off the device and unplug the power adaptor from the outlet.

## Before using the unit

### **ATTENTION**

- Please check the massage cloth cover for damage or tear before use. Do not use it if the cloth is torn. Please return the product to the after-sales service for inspection and repair.
- Please confirm that the operation buttons and timers are working properly before use.
- If you have not used for some time, please read the instruction manual carefully again, make sure the machine is normal, and then use after running safely.

#### **Inner cover of massager**

Please check if the cloth cover is damaged.


# Product Specifications

Model	FeelYou2 (Series JA)
Name	Massage Cushion
Power input	DC12V --- 2.5A
Power consumption	30 W
Default Time	20 mins


Short time operation 20 minutes, resting time 30 minutes

**NOTE:** Certified manufacturer's adapter must be equipped in the use of this products.

## Massage Position

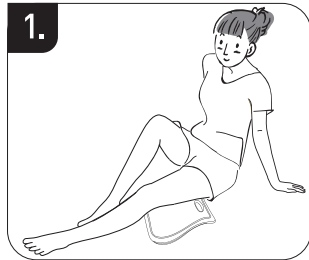Thighs

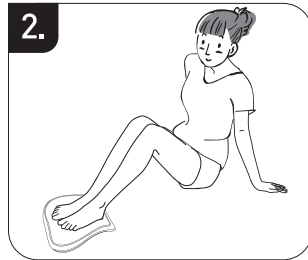
Feet


Abdomen

Back

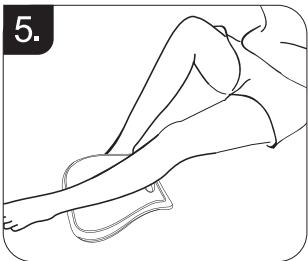calves

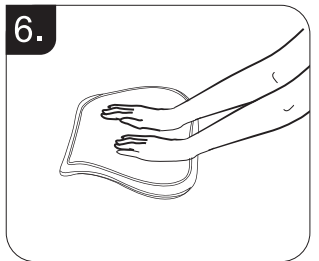hands

# How to clean

## 1. Cleaning

Before you clean, please make sure the power is switched off, and the power plug has been pulled out from the socket. Please use dry cloth to clean the product.


## 2. Maintenance

Keep the product clean and far from the wet environment.

**Clean when you feel cleaning is necessary.**

### Cleaning the main unit

- Wipe away dirt with a cloth soaked in soapy water, wring out the cloth before wiping.
- Do not use alcohol, paint thinner, benzine, or nail polish remover. Using these may damage the appliance or cause cracks or discoloration in the components.
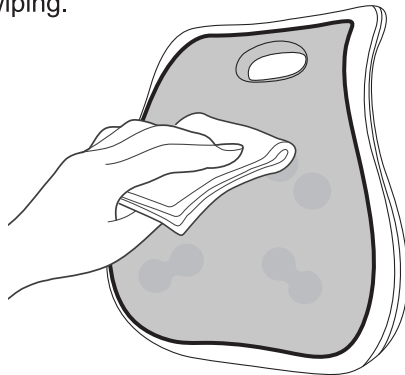
# Q&A

<b>Problem</b>	<b>Cause</b>	<b>Remedy</b>
Unable to start the device.	The cord is not connected properly.	Please push it harder until you heard the sound of "click".
The appliance stops during use.	The appliance automatically stops after approximately 20 minutes.	This does not indicate a problem.
The main unit does not heat up.	The appliance may not work beyond the recommended temperature.	Use under the recommended temperature of 15 °C - 35 °C

**If you still have any problem after using the methods above, please contact your local retailer.**

# Technological parameter

Product: Massage Cushion  
Model: FeelYou2 (Серия JA)  
Input: DC12V  $\overline{=}$  2.5A  
Power: 30W

Manufacturer: XIAMEN EMOKA HEALTHSCIENCE &  
TECHNOLOGY CO. LTD.

Production date: on the massager

Made in China

Manufacturer's address: D Building No. 100, JINFU ROAD,  
Chengnan Industrial Zone, TONG'AN District, Xiamen, China  
Phone: 0086-755-25122982

## **General importer in Russia:**

OOO "Max Concept", OGRN 1147453008532  
Moscow, Lyubertsy district, Kotelniki, 1st Pokrovsky  
PR-d, 5  
phone 8-800-700-14-29

[www.ergonova.ru](http://www.ergonova.ru)  
[info@ergonova.ru](mailto:info@ergonova.ru)

Full list of representative offices: <http://ergonova.ru/about/stores/> all copyrights to texts and images owned by Ergonova GmbH and may not be used otherwise only with the consent of the copyright holder.

Copyright © 2012.

Ergonova GmbH

All rights reserved.